

Penlon Prima 460 Anaesthesia System

ANAESTHESIA SOLUTIONS

- ♦ Electronic flow display and conventional flowmeter options
- ♦ Fully integrated ventilator and absorber
- ♦ Designed and manufactured in the UK
- ♦ Selectatec® compatible backbar with interlock
- ♦ HIS connectivity
- ♦ Intuitive design
- ♦ Low life costs



Penlon ♦

Penlon Prima 460 Anaesthesia System



All the features and options you need to configure a system to your exact specification

- 1 Electronic flow display and conventional flowmeter options are available
- 2 Six ventilation modes
- 3 New top shelf with secure GCX™ mounting system for patient monitors
- 4 SP M5 Patient Monitor
- 5 Territory-specific electrical outlet options
- 6 Selectatec® compatible backbar (two or three stations), with interlock on three station option
- 7 Up to four cylinders
- 8 Illuminated work space with pull-out writing surface
- 9 GCX™ compatible aluminium uprights for additional accessory mounting
- 10 Large capacity drawer units
- 11 Integrated CO₂ absorber and bellows unit with ventilator interface

The Penlon Prima 460 is a fully featured anaesthesia system providing the ideal solution for today's busy operating room

Clinician-focused choices and benefits, including electronic flow display options, generous work space, a functional top shelf module and GCX compatible monitor mount.

Platform

The Prima 460 Anaesthesia System has evolved from a strong core specification, including integrated AV-S Ventilator and A200SP Absorber.

- Fully compliant to ISO 80601-2-13 Standard and to the Restrictions of Hazardous Substances (RoHS) Directive.

Flow control and visualisation

The Prima 460 offers a wide range of options:

- Standard/cascade flow tubes
- Cascade flow tubes with electronic flow display providing numeric representation of gas flow

All models utilise a proven conventional control system with a mechanical anti-hypoxic system giving enhanced patient safety.

Models with an electronic flow display also have Hospital Information Systems (HIS) capability*.

Options and accessories

Wide choice of territory-specific electrical power outlets, and a forward facing socket. Penlon also supplies an AGS system (anaesthetic gas scavenging), a patient cable management arm, an oxygen therapy flowmeter, and a side mounted suction controller kit.

Patient monitoring

12.1" colour display with up to 12 waveforms, available with touchscreen control; up to 13 module slots (optional satellite module rack required); Nellcor™ or Masimo™ SpO₂ module options; Multi-parameter modules; Multi-Gas module with optional O₂ and BIS; and Sidestream, Mainstream or Microstream CO₂ modules.

System components

A reliable platform, combining advanced features and value for money. The Prima 460 is easy to use and maintain, with proven performance.

1 AV-S Ventilator

An easy to use, multifunction anaesthesia ventilator, designed for adult, paediatric and neonatal patient profiles. Standard specification includes three advanced spontaneous support modes and electronic PEEP. The AV-S has integrated oxygen monitoring and spirometry, plus seamless integration with both the absorber and anaesthetic machine, and HIS capability*.

2 A200SP Absorber

A high performance absorber with a ventilator interface as standard that provides ventilator mode switching, triggered by the bag/ventilator control. The unit has a quick-release canister, and the main components are autoclavable. A heater unit is offered as an option.

3 Sigma Delta Vaporizer

Award winning, high performance market leader, offering multi-choice fill system options.

4 SP M5 Patient Monitor


A lightweight and compact patient monitor with a 12.1" colour display, extensive parameter options, rechargeable Li-ion battery and optional wireless networking.

Maintenance and after-sales support

Penlon is committed to a successful, long term relationship with all our customers. Comprehensive warranty provides user peace of mind and after-sales support.

Additional services and warranties can be purchased to meet your particular needs.

*Relevant system interface is required

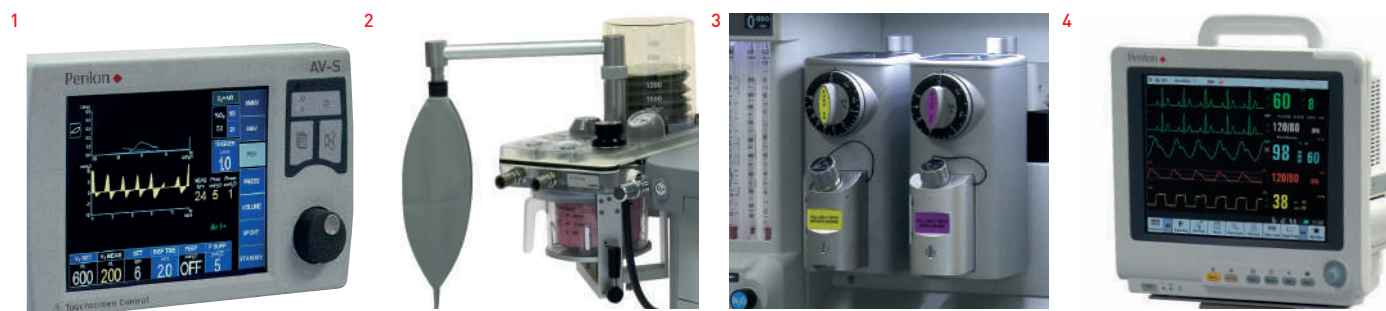


About Penlon ♦

Penlon was founded in 1943 by personnel from the Department of Anaesthesia at Oxford University. One of the first products was the Macintosh Laryngoscope, then a revolutionary design, and still the most widely used today, invented by the late Sir Robert Macintosh, Professor of Anaesthetics.

Today Penlon continues to design, engineer and build high quality anaesthesia products at its UK operations headquarters. The company is proud to have over 70 years' dedicated experience, many awards for product design, and an impressive four Queen's Awards for Enterprise, one for 'Innovation' and three for 'International Trade'.

Penlon devices feature intuitive user interfaces that require minimal operator training, putting clinicians in control, enabling them to focus on what is most important – patient safety and wellbeing.

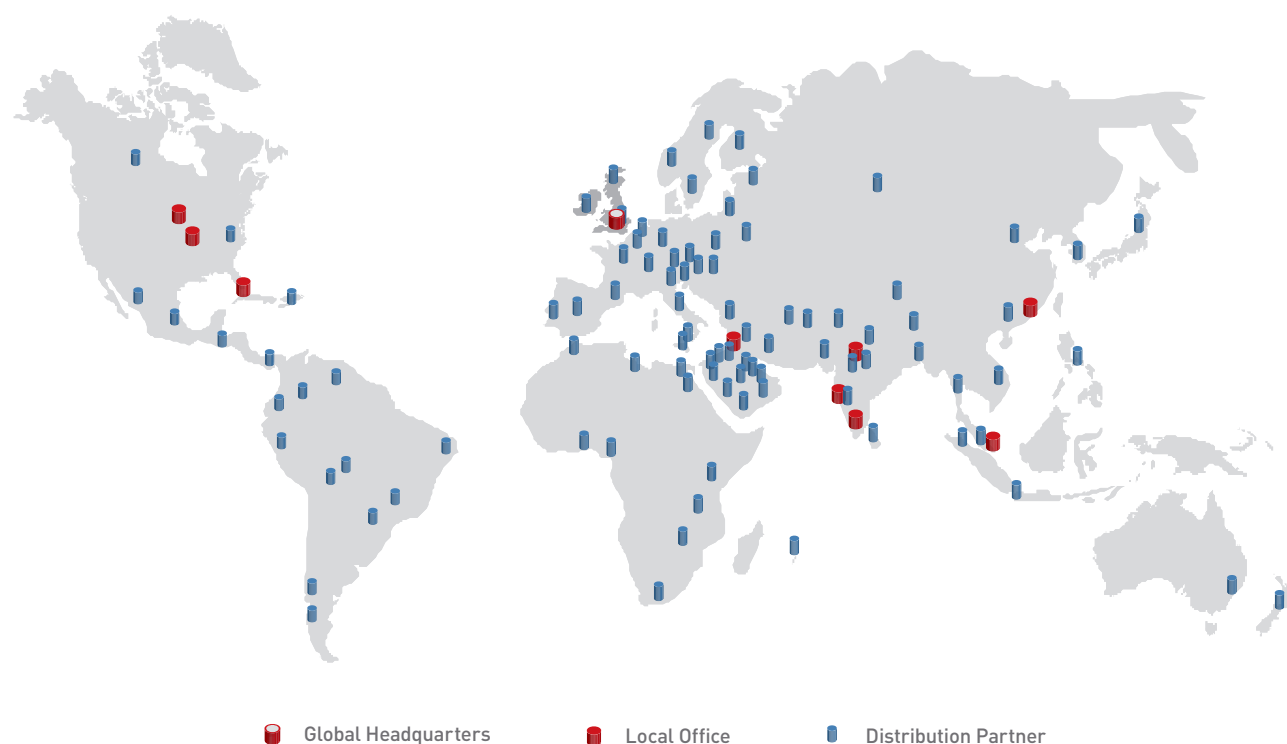

Service Helpline

1800-425-2355

CERTIFIED ISO 13485 : 2016 COMPANY

BPL Medical Technologies Private Limited

11th KM, Bannerghatta Road, Arakere, Bangalore 560 076

Toll Free: 1800-4252355

Website: www.bplmedicaltechnologies.com

For Enquiries: sales.medical@bpl.in

CIN: U33110KA2012PTC067282

Follow Us on



© 2018 BPL Medical Technologies Private Limited. All rights reserved. BPL Medical Technologies Private Limited reserves the right to make changes in product features, specifications, aesthetics and/ or to discontinue the same at any time without notice or obligation